

NURSE CALL ACCESSORIES AND FALLS PREVENTION

Need Help?

For assistance to choose the right product for your requirements do not hesitate to call one of our friendly sales team on 03 9588 0392 or email sales@safelife.com.au

Our falls prevention and nurse call accessories are compatible with all brand nurse call systems:

Smart Caller, Austco, Questek, Advanced Care, Xacom, Care Tech, Sedco, Smart Link & more.

Safe Life ABN: 67 637 195 941 Website: www.safelife.com.au

CONTENTS

FALLS PREVENTION HARDWIRED KITS	4 - 6
FALLS PREVENTION HARDWIRED REPLACEMENT PARTS	7 - 11
"OUT OF BED" SENSOR BEAMS	12
FALLS PREVENTION WIRELESS KITS	13 - 14
FALLS PREVENTION WIRELESS REPLACEMENT PARTS	15 - 18
ADAPTORS AND CORDS	19 - 20
NURSE CALL ACCESSORY PRODUCTS, PENDANTS	
AND BELL PRESS	21 - 23
PAGERS & PAGING KITS	24 - 25
WANDERING ALERT	26 - 27
INCONTINENCE & SEIZURE DETECTION	27

Please note:

Images and information in this catalogue were true and correct at the time of printing. Images, product specifications and price are subject to change without notice.

For detailed terms and conditions of sale, warranty and goods return arrangements - please contact Safe-Life at Alerting Devices Australia Pty Ltd, or your local agent/distributor.

FALLS PREVENTION HARDWIRED KITS



A-SL-BPHW-FOLDABLE-KIT

SAFE LIFE FOLDABLE HARDWIRED BED **EXIT KIT - FACILITY USE**

The Safe Life Bed Exit Kit comes complete with all parts for immediate Nurse Call use. This particular bed pad lays under the mattress protector/Kylie. When the resident lifts themself off the bed pad, the carers will be alerted via the nurse call system.

KIT CONTAINS:

- 1 x Safe Life Foldable Bed Exit Pad (A-SL-BPHW-01)
- 1 x Safe Life Hardwired monitor (A-SL-MON-HW)
- 1 x Power supply (A-CBM/BME-PP)
- 1 x Male to male cord (A-GCORD-2M)
- 1 x Double adaptor (A-ADAP-02-SM)
- 1 x Monitor soft cover and brackets





A-SL-CPHW-KIT-FAC

SAFE LIFE HARDWIRED CHAIR EXIT **KIT - FACILITY USE**

The Safe Life Chair Exit Kit comes complete with all parts for immediate Nurse Call use. Place the pad on the chair or seat. When the resident lifts themselves off the chair pad, the carers will be alerted via the nurse call system.

KIT CONTAINS:

- 1 x Safe Life Chair Exit Pad (A-SL-CPHW-01)
- 1 x Safe Life Hardwired monitor (A-SL-MON-HW)
- 1 x Power supply (A-CBM/BME-PP)
- 1 x Male to male cord (A-GCORD-2M)
- 1 x Double adaptor (A-ADAP-02-SM)
- 1 x Monitor soft cover and brackets



A-BPSC-KIT

CURA1 HARDWIRED CHAIR EXIT KIT - FACILITY USE

The Cura1 Chair Exit Kit comes complete with all parts for immediate Nurse Call use. Place the pad on the chair or seat. When the resident lifts themself off the chair pad, the carers will be alerted via the nurse call system.

- 1 x Cura1 Hardwired Chair Exit Pad (A-CPHW-01)
- 1 x Cura1 Hardwired monitor (A-BME-03)
- 1 x Power supply (A-CBM/BME-PP) 1 x Male to male cord (A-GCORD-2M)
- 1 x Double adaptor (A-ADAP-02-SM)
- 1 x Monitor soft cover and brackets

FALLS PREVENTION HARDWIRED KITS



A-BPLB-KIT-FACILTY

CURA1 HARDWIRED BED EXIT KIT - FACILITY USE

The Cura1 Bed Exit Kit comes complete with all parts for immediate Nurse Call use. This particular bed pad lays under the mattress protector/Kylie. When the resident lifts themself off the bed pad, the carers will be alerted via the nurse call system.

KIT CONTAINS:


- 1 x Cura1 Hardwired Bed Exit Pad (A-BPLB-01)
- 1 x Cura1 Hardwired monitor (A-BME-03)
- 1 x Power supply (A-CBM/BME-PP)
- 1 x Male to male cord (A-GCORD-2M)
- 1 x Double adaptor (A-ADAP-02-SM)
- 1 x Monitor soft cover and brackets



CURA1 HARDWIRED UNDER MATTRESS BED EXIT SENSOR KIT - FACILITY USE

The Cura1 Bed Exit Kit comes complete with all parts for immediate Nurse Call use. This particular bed pad lays under the mattress itself. When the resident lifts themself off the bed pad, the carers will be alerted via the nurse call system.

- 1 x Cura1 Hardwired Under mattress Bed Exit Pad (A-HW-MATTRESS-PAD)
- 1 x Cura1 Hardwired monitor (A-BME-03)
- 1 x Power supply (A-CBM/BME-PP)
- 1 x Male to male cord (A-GCORD-2M)
- 1 x Double adaptor (A-ADAP-02-SM)
- 1 x Monitor soft cover and brackets

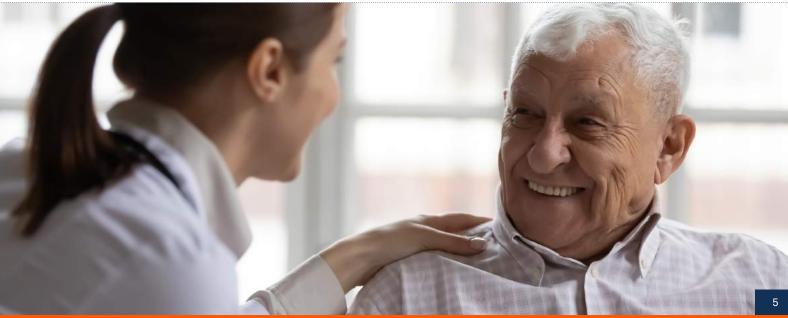

FALLS PREVENTION HARDWIRED KITS





A-SL-BP-CON-KIT

SAFE LIFE IN-LINE MONITOR BED EXIT KIT - FACILITY USE

Plugs direct into nurse call system. The Safe Life In-Line Monitor Bed Exit Kit comes complete with all parts for immediate Nurse Call use. This bed pad lays under the mattress protector/Kylie. No bulky monitor, just a neat, small in-line monitor. When the resident lifts themself off the bed pad, the carers will be alerted via the nurse call system. 2 second delay = less false alarms.

KIT CONTAINS:

- 1 x In-Line Monitor (A-SL-CON-M)
- 1 x In-Line Monitor Bed Exit Pad (A-SL-BP-CON-01)





A-SL-CP-CON-KIT

SAFE LIFE IN-LINE MONITOR CHAIR EXIT KIT - FACILITY USE

Plugs direct into nurse call system. The Safe Life In-Line Monitor Chair Exit Kit comes complete with all parts for immediate Nurse Call use. This chair pad is to be placed on the chair/seat. No bulky monitor, just a neat, small in-line monitor. When the resident lifts themself off the chair pad, the carers will be alerted via the nurse call system. 2 second delay = less false alarms.

KIT CONTAINS:

- 1 x In-Line Monitor (A-SL-CON-M)
- 1 x In-Line Monitor Chair Exit Pad (A-SL-CP-CON-01)





SAFE LIFE IN-LINE MONITOR UNDER MATTRESS BED EXIT SENSOR KIT - FACILITY USE

Plugs direct into nurse call system. The Safe Life In-Line Monitor U/mattress Bed Exit Kit comes complete with all parts for immediate Nurse Call use. This bed pad lays under the mattress itself. No bulky monitor, just a neat, small in-line monitor. When the resident lifts themself off the mattress, the carers will be alerted via the nurse call system. 2 second delay = less false alarms.

- 1 x In-Line Monitor (A-SL-CON-M)
- 1 x In-Line Monitor Chair Exit Pad (A-SL-CP-CON-01)

FALLS PREVENTION HARDWIRED REPLACEMENT PARTS







A-MAXIMAT-FOLDABLE A-MAXIMAT-FOLDABLE-BEIGE

SAFE LIFE FOLDABLE FLOOR SENSOR MAT HARDWIRED

The foldable floor pressure mat hard-wired is designed to immediately alert that a person has exited the bed.

- This mat can be folded into 3 for easy storage and shipping
- Durable quality vinyl
- · Foldable for storage and shipping
- Plugs straight into nurse call point via 6.5mm jack plug breakcord
- Comes with break cord and built in double adapter + spares available Ex Stock
- Measures unfolded 119 cm x 60 cm
- Folded size is 60 cm x 40 cm

Compatibility options available to suit all nurse call connection configurations – Please specify call point connection type on orders.







A-MAXIMAT-01

CURA1 FLOOR SENSOR MAT HARDWIRED

This large hardwired floor pressure mat is designed to immediately alert that a person has exited the bed.

- Large (1200 x 600mm) pressure mat
- · Anti-slip 'Yoga Top'
- Waterproof and washable using a low profile, non-slip material
- Breakcord with jack plug
- A double adapter (See section III) may be required if a bedside pendant cord is also in use.

Compatibility options available to suit all nurse call connection configurations – Please specify call point connection type on orders.



A-SL-MAXI MAT-02

SAFE LIFE FLOOR SENSOR MAT HARDWIRED

Much like the Cura1 hardwired floor mat, this large Safe Life model floor pressure mat is hard-wired and designed to immediately alert that a person has exited the bed.

- Large (1200 x 600mm) pressure mat
- Anti-slip 'Yoga Top'
- Waterproof and washable using a low profile, non-slip material
- Connects into Nurse Call Point via breakcord with 6.5mm jack plug and built in double adaptor

Compatibility options available to suit all nurse call connection configurations – Please specify call point connection type on orders

FALLS PREVENTION HARDWIRED REPLACEMENT PARTS



A-STANDARD-MAT-01

SAFE LIFE STANDARD FLOOR SENSOR **MAT - HARDWIRED**

This smooth vinyl floor sensor mat is one of our most popular.

- Measures 120 x 61 cm
- Waterproof and Washable
- Comes with break cord and built in double adapter + spares available Ex Stock
- 12 Month warranty

Compatibility options available to suit all nurse call connection configurations - Please specify call point connection type on orders.



A-CRASHMAT-CURA



CURA1 CRASH MAT WITH SENSOR - FACILITY USE

This floor crash mat is designed for residents who tend to fall out of the bed. Once impact with the mat is detected, this mat not only breaks the fall but will send an alarm through the nurse call system.

- · Connects to nurse call point
- Open dimensions 1800 x 750mm
- Folded dimensions 900 x 750mm
- 1 Year life expectancy
- Recommended weight range for reliable operation: 30kg to 200kg

Compatibility options available to suit all nurse call connection configurations - Please specify call point connection type on orders.



CRASH MAT (NO NURSE CALL SENSOR) - 50MM X 600MM

Place the Safety Crash Mat on the floor at the residents bed-side. Should the resident fall out of the bed the safety mat will help stop injuries by softening the impact.

- Foldable: made from three sections of safety foam that are then covered by hinged nylon
- The underside of the mat features a unique textured surface that prevents that mat from easily slipping if someone stands on it by mistake.
- Includes a reflective strip that runs along the edge of the mat. This strip will glow in the dark or even when the lights are low to make



A-SL-BPHW-01

SAFE LIFE FOLDABLE BED PAD - HARDWIRED

This hardwired bed pad is foldable for easy storage and shipping.

- Measuring 76cm x 52cm unfolded
- · Measuring 52cm x 38cm when folded
- The A-SL-HMON bedside control unit/monitor is a necessary accessory
- Placed under the fitted sheet or mattress protector at shoulder position
- Built in breakcord with RJ connection for monitor connectivity
- · Works with the Safe Life & Cura1 Monitors



A-BPLB-01

CURA1 BED PAD - HARDWIRED

This pressure sensitive bed exit sensor is usually placed under the Kylie mattress protector at shoulder position.

- Measuring 760 mm x 510 mm
- The A-BME-03 bedside control unit is a necessary accessory
- Comes with USB to RJ breakcord for monitor connectivity
- Works with the Safe Life & Cura1 Monitors



A-HW-MATPAD

CURA1 UNDER MATTRESS BED PAD - HARDWIRED

This hardwired bed exit pad is to be placed underneath the mattress rather than under the mattress protector.

- · Dimensions 760 x 510mm.
- Recommended weight range for reliable operation: 40kg to 200kg
- The A-BME-03 bedside control unit is a necessary accessory.
- · Comes with USB to RJ breakcord for monitor connectivity
- · Works with the Safe Life & Cura1 Monitors



A-SL-CPHW-01

SAFE LIFE CHAIR EXIT PAD - HARDWIRED

The Safe Life chair pad offers the same falls detection functionality as the Safe Life bed pad in a smaller chair friendly size. It will trigger an alarm along with audible tone in the event of the resident leaving or falling from the chair.

- Measures 38 x 19 cm
- · Low Cost and long lasting
- · Works with the Safe Life & Cura1 Monitors
- Considered a consumable item

FALLS PREVENTION HARDWIRED REPLACEMENT PARTS





A-SL-BP-CON-01

SAFE LIFE IN-LINE MONITOR BED PAD - FACILITY USE (IN-LINE MONITOR NOT INCLUDED)

The latest release of the Safe Life range is the In-Line Monitor Bed pad.

This Bed pad does not require a standard monitor for operation; however, it does need an in-line monitor (A-SL-CON). This makes it extremely cost effective. Simply plug directly into the bedside call point/button. This is a foldable bed bad, as seen folded to the left.



CURA1 CHAIR EXIT PAD - HARDWIRED

Where a chair is in close proximity to a call-point this chair exit sensor can be hard-wired to its BME-03 controller (listed below). It will trigger an alarm along with audible tone in the event of the resident leaving or falling from the chair.

- Measures 35 x 38 cm
- The A-BME-03 bedside control unit/monitor is a necessary accessory
- · Works with the Safe Life & Cura1 Monitors



SAFE LIFE IN-LINE MONITOR UNDER MATTRESS BED PAD - FACILITY USE (IN-LINE MONITOR NOT INCLUDED)

The latest release of the Safe Life range is the In-Line Monitor Under Mattress Bed pad.

This Bed pad does not require a standard monitor for operation; however, it does need an in-line monitor (A-SL-CON). This makes it extremely cost effective. Simply plug directly into the bedside call point/button.



SAFE LIFE IN-LINE MONITOR CHAIR PAD - FACILITY USE (IN-LINE MONITOR NOT INCLUDED)

The latest release of the Safe Life range is the In-Line Monitor chair pad.

This chair pad does not require a standard monitor for operation; however, it does need an in-line monitor (A-SL-CON). This makes it extremely cost effective. Simply plug directly into the bedside call point/button.





A-SL-MON-HW

SAFE LIFE BED-EXIT, CHAIR-EXIT AND FLOOR MAT MONITOR - HARDWIRED

This low-cost monitor is an optional accessory for floor sensor mat use but is essential for a hardwired bed and chair exit sensors

The A-SL-HMON can be powered by a 9v battery (included) or by power pack.

- Default reset Which requires the two buttons on the side of the monitor to be pressed simultaneously in order to avoid residents attempting to reset the alarm
- Volume control
- · Adjustable time delay limiting false alarm calls
- Includes a socket for connection to the bedside call-point
- Temporary disarm which allows for a resident to be moved from a bed or chair by the carer without activating an alarm or making a noise.
- For operation through the nurse call system a A-GCORD is required





CURA1 BED-EXIT AND CHAIR EXIT MONITOR - HARDWIRED

This monitor is an optional accessory for hardwired floor sensor mat use but is essential for a hardwired bed and chair exit sensors.

The A-BME-03 Hardwired monitor can be powered by a battery (included) or by power pack.

- Reset button
- · Alarm delay adjustment
- Volume control for the audible alert Pause mode to attend to a patient without triggering the alarm
- Includes a socket for connection to the bedside call-point
- For operat

A-BME-03





SAFE LIFE IN-LINE MONITOR - FACILITY USE

This In-Line Monitor is for our new range of hardwired undermattress, bed and chair exit pads. Simply plug one of our new bed/chair pads into the call point via the In-Line Monitor lead. This "plug and play" design makes for a more affordable solution.

2 second delay = less false alarms.

- · 18 Months battery life expectancy
- User replaceable batteries
- Low Battery warning

Compatibility options available to suit all nurse call connection configurations.

PIR SENSOR BEAMS



OUT OF BED CURTAIN SENSOR BEAM - SINGLE PIR

This small Curtain PIR fits to the wall at the head of the bed in a purpose-built wall mount bracket. It is easy to install and can be moved from room to room with ease.

It projects an invisible beam to the end of the bed. When the resident's legs break this beam an alarm call is generated. Connects to all nurse call systems. Comes complete with wall mount bracket.


All cables are at the head of the lead eliminating trip hazards.

Compatibility options available to suit all nurse call connection configurations.

For a video demonstration please visit us at https://safelife.com.au/product/single-bed-exit-pir/





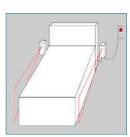





OUT OF BED CURTAIN SENSOR BEAM - DUAL PIR

This PIR has dual interconnected sensors. One PIR sensor each side of the bed makes it suitable for beds that can be accessed from either side (not against a wall).

Comes with a power pack wall mounting brackets and an in-Line ON/OFF switch located where it connected to the call point. Comes complete with 2 wall mount brackets.

Compatibility options available to suit all nurse call connection configurations.

For a video demonstration please visit us at https://safelife.com.au/product/single-bed-exit-pir/



Cura1 Wireless PIR Sensor Beam Kit with Multiport Receiver – Facility Use

Easily mountable and simple to install, the wireless PIR curtain sensor beam is the perfect addition to all set-ups where the resident is a falls risk or may need to be monitored. The multiport receiver will plug into the bedside call point and send an alarm through to the nurse call system when the iBeam

PIR senses that someone has walked through a door way, or is attempting to get out of bed.

The Cordless Multiport receiver can monitor up to 5 wireless

devices. (2 x pads and 3 x other transmitters eg. PET)



Cura1 Wireless PIR Sensor Beam Kit – Facility Use

Easily mountable and simple to install, the wireless PIR curtain sensor beam is the perfect addition to all set-ups where the resident is a falls risk. The Cordless Falls Monitor will plug into the bedside call point and send an alarm through to the nurse call system when the iBeam PIR senses that someone has walked through a door way, or is attempting to get out of bed.

Unlike the Multiport Receiver kit above, an interconnecting cord from the Cordless Falls Monitor into the bedside call point means this kit can be used as a standalone system or integrated into the nurse call/call bell system. The Cordless Falls Monitor can monitor up to 6 wireless Cura1 devices — which means you can add 5 more sensors to this kit (chair exit sensor pad, bed exit sensor pad, wireless help button and more).

FALLS PREVENTION WIRELESS KITS





A-CL-MULTIPORT-FM-KIT

CURA1 WIRELESS FLOOR SENSOR MAT KIT - WITH MULTIPORT RECEIVER - FACILITY USE

This cordless floor sensor mat kit is designed to alert that a person has exited a bed/chair or entered a room via the multiport receiver cordless monitoring device. The multiport receiver plugs directly into the nurse call system call point.

Connection is Tip & Sleave. Tip & Ring available on request.

KIT CONTAINS:

- 1 x Cura1 Wireless Transmitter (A-CL433-TX-5V)
- 1 x Cura1 Wireless Floor Sensor Mat Pad (A-MAXIMAT-CL)
- 1 x Cura1 Wireless Multiport receiver (Mini cordless monitor plugs direct to nurse call) (A-CL-MULTIPORT)



A-CL-MULTIPORT-BP-KIT

CURA1 WIRELESS BED EXIT KIT - WITH MULTIPORT RECEIVER - FACILITY USE

This cordless bed exit kit is designed to alert that a person has exited the bed via the multiport receiver cordless monitoring device. The multiport receiver plugs directly into the nurse call system call point.

Connection is Tip & Sleave. Tip & Ring available on request.

KIT CONTAINS:

- 1 x Cura1 Wireless Transmitter (A-CL433-TX-5V)
- 1 x Cura1 Wireless Bed Pad (A-CBP-01)
- 1 x Cura1 Wireless Multiport receiver (Mini cordless monitor plugs direct to nurse call) (A-CL-MULTIPORT)



CURA1 WIRELESS UNDER MATTRESS BED EXIT - WITH MULTIPORT RECEIVER - FACILITY USE

This cordless under mattress bed exit kit is designed to alert that a person has exited the bed via the multiport receiver cordless monitoring device. The multiport receiver plugs directly into the nurse call system call point. Connection is Tip & Sleave. Tip & Ring available on request.

KIT CONTAINS:

- 1 x Cura1 Wireless Transmitter (A-CL433-TX-5V)
- 1 x Cura1 Wireless Under Mattress Sensor Pad (A-CL-MATTRESS-PAD)
- 1 x Cura1 Wireless Multiport receiver (Mini cordless monitor plugs direct to nurse call) (A-CL-MULTIPORT)



CURA1 WIRELESS CHAIR EXIT KIT - WITH MULTIPORT RECEIVER - FACILITY USE

This cordless chair exit kit is designed to alert that a person has exited the chair or seat via the multiport receiver cordless monitoring device. The multiport receiver plugs directly into the nurse call system call point. Connection is Tip & Sleave. Tip & Ring available on request.

- 1 x Cura1 Wireless Transmitter (A-CL433-TX-5V)
- 1 x Cura1 Wireless Chair Pad (A-CCP-01)
 - 1 x Cura1 Wireless Multiport receiver (Mini cordless monitor plugs direct to nurse call) (A-CL-MULTIPORT)

FALLS PREVENTION WIRELESS KITS



A-MAXI CL KIT



CURA1 WIRELESS FLOOR SENSOR MAT KIT - FACILITY USE

This large cordless floor sensor mat kit is designed to alert that a person has exited the bed/seat or exited a room via the cordless falls monitor. The cordless falls monitor plugs into the nurse call system call point via the male to male cord supplied. The monitor can also be used as a portable standalone system.

KIT CONTAINS:

- 1 x Cura1 Wireless Transmitter (A-CL433-TX-5V)
- 1 x Cura1 Wireless Floor Sensor Mat (A-MAXIMAT-CL)
- 1 x Cura1 Wireless Monitor (A-CBM-02)
- 1 x Power supply (A-CBM/BME-PP)
- 1 x Male to male cord (A-GCORD-2M)
- 1 x Double adapter (A-ADAP-02-SM)
- 1 x Monitor soft cover and brackets



A-CBP-KIT-FACILITY



CURA1 WIRELESS BED EXIT SENSOR KIT - FACILITY USE

This cordless bed exit kit is designed to alert that a person has exited the bed via the cordless falls monitor. The cordless falls monitor plugs into the nurse call system call point via the male to male cord supplied. The monitor can also be used as a portable standalone system.

KIT CONTAINS:

- 1 x Cura1 Wireless Bed Exit Pad (A-CBP-01)
- 1 x Cura1 Wireless Transmitter (A-CL433-TX-5V)
- 1 x Cura1 Wireless Monitor (A-CBM-02)
- 1 x Power supply (A-CBM/BME-PP)
- 1 x Male to male cord (A-GCORD-2M)
- 1 x Double adapter (A-ADAP-02-SM)
- 1 x Monitor soft cover and brackets



A-CCP-KIT

CURA1 WIRELESS CHAIR EXIT SENSOR KIT - FACILITY USE

This cordless chair exit kit is designed to alert that a person has exited the chair via the cordless falls monitor. The cordless falls monitor plugs into the nurse call system call point via the male to male cord supplied. The monitor can also be used as a portable standalone system.

KIT CONTAINS:

- 1 x Cura1 Wireless Chair Exit Pad (A-CCP-01)
- 1 x Cura1 Wireless Transmitter (A-CL433-TX-5V)
- 1 x Cura1 Wireless Monitor (A-CBM-02)
- 1 x Power supply (A-CBM/BME-PP)
- 1 x Male to male cord (A-GCORD-2M)
- 1 x Double adapter (A-ADAP-02-SM)
- 1 x Monitor soft cover and brackets



CURA1 WIRELESS UNDER MATTRESS BED EXIT SENSOR KIT - FACILITY USE

This cordless bed exit kit is designed to alert that a person has exited the bed via the cordless falls monitor. This Particular pad is positioned under the mattress, rather than just under the mattress protector/kylie. The cordless falls monitor plugs into the nurse call system call point via the male to male cord supplied. The monitor can also be used as a portable standalone system.

- 1 x Cura1 Cordless Under mattress Bed Exit Pad (A-CL-MATTRESS-PAD)
- 1 x Cura1 Cordless Transmitter (A-CL433-TX-5V)
- 1 x Cura1 Cordless Monitor (A-CBM-02)
- 1 x Power supply (A-CBM/BME-PP)
- 1 x Male to male cord (A-GCORD-2M)
- 1 x Double adapter (A-ADAP-02-SM)
- 1 x Monitor soft cover and brackets

FALLS PREVENTION WIRELESS REPLACEMENT PARTS





A-CBM-02

CURA1 RECEIVER/CONTROLLER FOR WIRELESS DEVICES

Multiple wireless device monitoring capability.

- Monitor up to 7 devices simultaneously a single Premium pad/mat and up to 6 wireless devices
- · Comes with durable ABS plastic enclosure and Caregiver Reset Key.
- Alarm memory and auto re-arm
- · Amber LED visual alarm indicator
- Caregiver secure reset option
- Battery operated or optional external plug pack, up to 6 months battery life (nominal)
- 3 selectable alarm tones and adjustable alarm volume
- "Pause" feature allows caregiver to attend to the patient without setting the alarm off
- Alerts for lost nurse call connection, lost AC power supply & low battery
- Dimensions 120 x 77 x 40mm
- Up to 30m range when used with Wireless bed pads and chair pads

Item requires GCORD to allow for nurse call operation





A-CL-MULTIPORT (Facility use only)

CURA1 MULTIPORT RECEIVER FOR WIRELESS DEVICES

Multiple wireless device monitoring capability.

Much Like the above device, the Multiport Receiver is a more compact, simplified version which will plug directly into a call point.

- Monitor up to 5 wireless devices and 1 hardwired device simultaneously (including pendant cord OR hardwired floor mat)
- · Cord removal alarm
- Dimensions 78 x 46 x 33mm
- · LED status & warning indicator lights
- Plugs directly into most nurse call systems Stereo Tip & Sleeve config, however on request may be modified and supplied Tip & Ring at an additional cost.
- Battery operated or optional external plug pack, 3 months battery life (nominal)
- · Up to 8-12m range when used with wireless devices





CURA1 MOUNTED MULTIPORT RECEIVER FOR WIRELESS DEVICES

- Installed on a bracket away from the nurse call panel so it is suitable for congested nurse call bed head panels where directly mounting on the panel isn't possible.
- Works with up to 5 wireless devices. (2 x pad & 3x other transmitters eg PET)
- Works with: Cordless Bed/Chair Pad or Floor Mat; Personal Emergency Transmitter (PET); iBeam; Wired Nurse Call Pendant; Wireless Magnetic Reed Switch; Wireless Waterproof Call Buttons.
- · LED Status and warning indicator lights.
- 3 Minute Pause Mode
- · Cord removal alarm.

FALLS PREVENTION WIRELESS REPLACEMENT PARTS



A-ACTIVECARE-PAGER

(Home Care)





- Large easy-to-read 7-digit LCD screen displays device type and ID.
- Can be programmed to monitor up to 9 wireless devices
- Up to 5 simultaneous alarms can be stored in the pager's "Call Review" memory
- Vibration and tone selection for hard of hearing.
- Out of range check.
- Wireless reset capability.
- Direct signalling from wireless devices for use in home care environments.
- Easily cleanable buttons



A-MAXIMAT-CL

CURA1 FLOOR SENSOR MAT - WIRELESS

This Wireless product is designed to immediately alert that a person has exited the bed.

- Large soft 'yoga' top mat
- Measures 120 cm x 60 cm
- Waterproof and washable using a low profile, non-slip material
- 'Out-rigger' press-stud pocket that will allow the insertion of an A-CL-433M-TX-V5 transmitter (necessary accessory)

Under this part number no transmitter is supplied. Use this part number to purchase replacement floor mats. Either the A-CBM-02, A-ACTIVECARE-PAGER, A-CL-MULTIPORT or A-CL-MULTIPORT-02 is a necessary monitoring accessory (not included)

Also see Complete nurse call kits





A-MAXI-MAT-433

CURA1 WIRELESS FLOOR SENSOR MAT INCL. TRANSMITTER

Same Floor mat sensor as above, however this part number includes the transmitter

Either the A-CBM-02, A-ACTIVECARE-PAGER, A-CL-MULTIPORT or A-CL-MULTIPORT-02 is a necessary monitoring accessory (not included) Also see Complete nurse call kits





CURA1 UNDER MATTRESS BED EXIT SENSOR PAD - WIRELESS

This Wireless bed exit pad is to be placed underneath the mattress rather than under the mattress protector.

- Pad dimensions: 800 x 290mm
- Needs A-CL433-TX-5V Transmitter (necessary accessory)
- Recommended weight range for reliable operation: 40kg to 200kg
- Material: Non allergenic PVC
- Under this part number no transmitter is supplied. Use this part number solely to purchase replacement pads.

Either the A-CBM-02, A-ACTIVECARE-PAGER, A-CL-MULTIPORT or A-CL-MULTIPORT-02 is a necessary monitoring accessory (not included) Also see Complete nurse call kitsAlso see Complete nurse call kits



CURA1 CHAIR EXIT SENSOR PAD - WIRELESS

This Wireless chair pad is used on chairs and wheel chairs of persons subject to falls risks. Can be used on paradise chairs, recliners and lift chairs.

- Dimensions 380 x 255mm
- To be used with A-CL433-TX-V5 Transmitter (necessary accessory)
- Recommended weight range for reliable operation: 30kg to 200kg
- Under this part number it includes no transmitter and is the part number to use solely for supply of replacement pads

Either the A-CBM-02, A-ACTIVECARE-PAGER, A-CL-MULTIPORT or A-CL-MULTIPORT-02 is a necessary monitoring accessory (not included)



A-CCP-01





CURA1 WIRELESS CHAIR PAD INCL. TRANSMITTER

Same chair exit pad as above, however this part number includes the transmitter

Either the A-CBM-02, A-ACTIVECARE-PAGER, A-CL-MULTIPORT or A-CL-MULTIPORT-02 is a necessary monitoring accessory (not included)



A-CBP-01



CURA1 BED EXIT SENSOR PAD - TRANSMITTER WIRELESS

This pressure sensitive bed exit sensor is usually placed under the Kylie mattress protector at shoulder position.

- Dimensions 760 x 510mm
- Needs A-CL433-TX-5V Transmitter (necessary accessory)
- Recommended weight range for reliable operation: 30kg to 200kg
- Under this part number no transmitter is supplied. Use this part number solely to purchase replacement pads.

Either the A-CBM-02, A-ACTIVECARE-PAGER, A-CL-MULTIPORT or A-CL-MULTIPORT-02 is a necessary monitoring accessory (not included)

Also see Complete nurse call kits



CURA1 WIRELESS BED PAD INCL

Same bed exit pad as above, however this part number includes the transmitter

Either the A-CBM-02, A-ACTIVECARE-PAGER, A-CL-MULTIPORT or A-CL-MULTIPORT-02 is a necessary monitoring accessory (not included)

FALLS PREVENTION WIRELESS REPLACEMENT PARTS



CURA1 WIRELESS TRANSMITTER

This Transmitter fits within the Wireless Floor, Bed and Chair sensor pads/mats transmitter enclosure and communicates with Cura1 cordless monitor, multiport receivers and

ActiveCare Pager

- It is a necessary accessory for Wireless pads/mats
- · 433MHz operating frequency
- Up to 50m range
- · 1-year battery life
- · Incorporates status/in-use LED indication
- · Incorporates on/off switch

Also see Complete nurse call kits



CURA1 PULL STRING MONITOR

- Used to detect a resident attempting to leave a chair or bed unassisted VIA a magnetic pull string and alligator clip
- When the resident pulls away from the direction of the monitor the monitor will alert with an audible tone
- Optional nurse call connectivity



SAFE LIFE PULL STRING MONITOR/SENSOR PAD MONITOR

- Used to detect a resident attempting to leave a chair or bed unassisted VIA a magnetic pull string and alligator clip
- When the resident pulls away from the direction of the monitor the monitor will alert with an audible tone
- Optional nurse call connectivity
- Can also connect and monitor 1 hardwired chair pad, bed pad or floor sensor mat
- · 3 selectable alarm tones



ADAPTORS AND CORDS



A-SL-BREAKCORD-01

SPARE BREAKCORD FOR SAFE LIFE HARDWIRED FLOOR MATS

- 6.35mm jack plug for nurse call wall plate connectivity
- Includes built in double adaptor

Compatibility options available to suit all nurse call connection configurations – Please specify call point connection type on orders



SPARE CURA1 BREAKCORD MONO

USB to stereo 6.35mm jack plug (MONO) for connectivity from Cura1 hardwired floor mat to nurse call wall plate





A-BREAKCORD-USB - T & R OR A-BREAKCORD-USB - T & S

SPARE CURA1 BREAKCORD – STEREO

USB to stereo 6.35mm jack plug (Tip and Ring) OR (Tip and Sleave) for connectivity from Cura1 hardwired floor mat to nurse call wall plate

Please specify T & S or T & R



A-BREAKCORD-USB-RJ10

SPARE CURA1 BREAKCORD – USB TO RJ10

USB to RJ10 plug for connectivity from Cura1 hardwired floor mats and Pads to Cura1 Monitor.





A-CBM & BME-PP

CURA1 POWER PACK FOR MONITORS – SUITS SAFE LIFE MONITOR ALSO

Plug pack for Cura1 and Safe Life falls monitors

ADAPTORS AND CORDS



TRIPLE ADAPTORS - LOW PROFILE 6.35MM - UNIVERSAL CONNECTION (T-T, R-R, S-S)

This range of triple adaptors connect to the call-point socket whilst also allowing connection of multiple devices. For example, the triple adaptor would accommodate the normal over-bed cord pendant plus perhaps a floor mat and a bed exit kit or any combination of 3 devices that may be required from time to time.



DOUBLE ADAPTORS - LOW PROFILE 6.35mm - Universal Connection (T-T, R-R, S-S)

This range of double adaptors connect to the call-point socket whilst also allowing connection of multiple devices. For example, the double adaptor would accommodate the normal over-bed cord pendant plus perhaps a floor mat or a combination of any two devices that may be required from time to time.

A-ADAP-02-S-1



PENDANT EXTENSION LEAD 1M

For use when extra cord length for pendants, floor mats, monitors etc. is needed.

Our standard connection is mono 6.35mm - If you require a different connection please advise our sales team



PENDANT EXTENSION LEAD 2-3M

For use when extra cord length for pendants, floor mats, monitors etc. is needed.

Our standard connection is mono 6.35mm - If you require a different connection please advise our sales team

NURSE CALL ACCESSORY PRODUCTS, PENDANTS AND BELL PRESS



INTERCONNECTING CORD WITH TWO 6.35mm MALE JACKS TO CONNECT CONTROLLER TO CALL-POINT SOCKET

Will interconnect a Safe Life or Cura1 controller/Monitor to the bedside call-point which uses a 6.35mm socket.

- · 2 Meters length
- · Low profile plugs to protect call point from knocks
- · 6.35mm plug on each end
- · Discreet white colour

Our standard connection is mono $6.35 \mathrm{mm}$ - If you require a different connection please advise our sales team



BARREL GRIP PENDANT 3 OR 5 METRES

This ergonomically designed pendant is sleek with a moulded grip handle for easy use amongst elderly residents. 3- and 5-meter length available – Please specify length.

6.35mm right angle plug.

The 3M barrel grip comes in both MONO OR STEREO – please specify connection type at time of purchase. 5M version is MONO only.





TEAR DROP HANDLE PENDANT CORD

3m pendant nurse call pendant. Stereo 6.35mm right angle plug- tip and ring.

A-PC731



A-SL-PRESSURE-PEND-TR A-SL-PRESSURE-PEND-TS A-SL-PRESSURE-PEND-MONO

PRESSURE PAD (PADDLE) PENDANT CORD

This pressure pad pendant cord (also known as the "Paddle Pendant") is for residents with limited mobility. Placing a light weight or pressure on the 80mm round disk will activate an alarm. Suitable for hand, arm and knee operation.

- · Soft ergonomic feel
- · Waterproof and dust proof IP67
- · Anti bacterial + Anti fungal properties for infection control
- · 6.3mm plug configurations to suit all nurse call system brands
- · Includes fastener clip
- · Standard 3 meter length
- · Braille for the visually impaired
- · Eyelet for mounting in awkward positions

NURSE CALL ACCESSORY PRODUCTS, PENDANTS AND BELL PRESS



A-SILICON-1-M-3M A-SILICON-1-T/S-3M A-SILICON-1-T/R-3M



SINGLE BUTTON WATERPROOF SILICON CORD PENDANT - GLOW IN THE DARK -

This sealed pendant is waterproof, suitable for dip sterilization, unbreakable and uses anti-bacterial and anti-fungal materials. Its ergonomic design suits left and right-hand grip and has a Braille imprint for visually impaired residents/patients.

Available in different connection configurations to suit all brand nurse call systems.



A-SILICON-2-TR/S A-SILICON-2-ST/SR



TWO BUTTON WATERPROOF SILICON CORD PENDANT - GLOW IN THE DARK -

Two button for call and light switch operation. Like all Safe Life pendants we have a version of this silicon pendant to suit all brands of nurse call systems. This sealed pendant is made with anti-bacterial and anti-fungal materials. It is waterproof and suitable for dip sterilization. Its ergonomic design suits left and right-hand grip and has a Braille imprint for visually impaired residents/ patients.

Available in different connection configurations to suit all brands of nurse call systems. Light switching function also determined by nurse call system capability.



A-BP-PET-304-GRA



BLUE-PHONE/HP3/4/5 WIRELESS PENDANT

- This is the standard issue waterproof pendant (as supplied with the Smart-Caller range of ILU Help Phones and Diallers).
- This pendant is also used to trigger Smart Caller IP and ISM Wireless call points remotely
- The pendant can be opened to change the battery (Panasonic lithium battery)



REMOTE WRIST BAND PENDANT

The remote wrist band transmitter will activate the Pendant Alert Kit (below).



PENDANT ALERT KIT - GIVE ANY ROOM WIRELESS PENDANT CAPABILITIES - FACILITY USE ONLY

The Remote Pendant Alert Kit allows for a resident to move freely within their own room and call for assistance with the wireless wrist style pendant. Simply connect the Pendant Alert to the resident's bedside call point and give the resident their wireless pendant.

When the resident presses the pendant the care staff will be notified through the facility's nurse call system.

Compatibility options available to suit all nurse call connection configurations

 Please specify call point connection type on orders. If not specified 6.35mm mono connection will be supplied.



CURA 1 WIRELESS DURESS NECK PENDANT

- · Large rubber button for easy use by elderly patients
- Standard Lanyard for patient to wear around neck
- · Works with single use hospital wristbands.
- · Changeable battery with projected 1-year life.
- Splash resistant

Either the A-CBM-02, A-ACTIVECARE-PAGER, A-CL-MULTIPORT or A-CL-MULTIPORT-02 is a necessary (Not Included)





CURA1 WIRELESS MAGNETIC REED SWITCH

- Sends wireless signal when door or window is opened
- Built-in call button extends usage options e.g. use it as a portable panic button
- Double sided tape included to for easy mounting
- · Battery included
- · LED alarm active indication

Either the A-CBM-02, A-ACTIVECARE-PAGER, A-CL-MULTIPORT or A-CL-MULTIPORT-02 is a necessary (Not Included)

PAGERS & PAGING KITS



LONG RANGE CALL POINT TO PAGER KIT

The Call Point to Pager kit is excellent for use in home care, medical centres, hospital wings, aged care facilities and other public building spaces. 100 meters range, and more with the use of a range extender, this system provides peace of mind for both carers, and those whom may call for assistance.

Each call point will have an individual ID message which transmits directly to the professional pager/s.

There is no limit to the number of call points or pagers that can be used. Also, extremely valuable for staff duress in public buildings such as libraries. The applications are endless.



LONG RANGE PENDANT TO PAGER KIT

The Pendant to Pager kit is excellent for use in home care, medical centres, hospital wings, aged care facilities and other public building spaces.

100 meters range, and more with the use of a range extender, this system provides peace of mind for both carers, and those whom may call for assistance. Each pendant will have an individual ID message which transmits directly to the professional pager/s.

There is no limit to the number of pendants or pagers that can be used. Also, extremely valuable for staff duress in public buildings such as libraries.





CURA1 HOME CARE PAGER KIT

The Home Care Pager Kit offers direct communication between resident and Carer's pager in the form of an alarm tone. This Cura1 pager kit is a portable mini system providing wireless button to pager alerting application.

It is packaged to include two wireless pendants each having a mounting bracket plus lanyard for flexible use.

The system is suitable for covering short distances (up to 30 meters) in the home.



A-PAGER-7950-VHF or UHF 2-LINE LARGE SCREEN PAGER

This model POCSAG alphanumeric belt pager offers a good quality robust pager. It includes holster, handbook and battery.

Having eight addresses (cap codes) it has 4 beep tone cadences to identify the incoming event and priority and includes flashing back-light and user selectable vibrator function.

A heavy-duty silicon holster and /or Bungee cord options are also available. Silicon Holsters also available

Silicon Holsters also available RRP \$30.25 INCL. GST

Please specify VHF or UHF when purchasing



UHF PENDANT

The miniature pendant is a high powered pendant perfect for the small SRS or nursing home, Medical Centre or independent living units.

The UHF one button pendant is operated by a changeable lithium battery. The call point can transmit a signal to a pager up to 100 meters (in a clear line of sight).

Range extenders are available which can increase transmission range.





ONE BUTTON WIRELESS UHF CALL POINT

The one button mobile call point is perfect for the small SRS or nursing home, Medical Centre or independent living units.

The UHF one button pendant is operated by a changeable lithium battery.

The call point can transmit a signal to a pager up to 100 meters (in a clear line of sight).

Pendant cords, floor mats, bed and chair pads and sensor beams can also be connected to detect of of bed or chair.

Range extenders are available which can increase transmission range.

A-SAL-GP-PC-01

WANDERING ALERT



GATEWATCH-01



THE SAFELIFE GATEWATCH SYSTEM

The Australian designed and manufactured GateWatch alerting system uses the International Bluetooth Technology for its Door Location 'Beacons' and it's wide range of Personal Tags.

This system not only accommodates safe containment of residents and patients that tend to wander but can also allow both local and/or central alerting of assets that are leaving the building.

For example, the local proximity beeping option can avoid the problem of staff accidentally carrying Pagers, Tablets, Mobile and DECT phones home.

FEATURES AND BENEFITS

- Suits internal or outside locations via its waterproof design
- Adjustable Tag distance proximity via internal switch settings
- Connection for external alert e.g. strobe light or remote buzzer
- Satellite (clone) units for multiple or complex door locations
- · Staff/Family override capability for authorized exit
- Contact closure to suit all brands of nurse call/Access Control systems
- Contacts provided for door latch/unlatch option
- Long battery Life Bluetooth Tags with replaceable battery
- Easy to install and designed for maintenance free operation





M-SAFELIFE-PET-4G-2

NEW GENERAGTION 4G SAFE LIFE PERSONAL ALARM - FOR USE IN THE COMMUNITY AND AT HOME

This small 'go anywhere' personal alarm is perfect for dementia sufferers, people with disabilities, children, lone workers and those who live alone.

This device works in the community as well as at home. Press the SOS button to alert your carers network or family when you need help. Includes falls detection and GPS location.

- Personalised out-going emergency SMS message and emergency phone calls.
- Free shipping
- Improved GPS positioning
- Improved Falls Detection
- Two-way voice contact
- SOS SMS messages with map location
- Smart docking station/ charging base included (includes BT pinpoint location for in the home)
- Monitored by 1-10 family/ friends' mobile phones
- You may apply 000 as your last emergency contact
- Geo Fence
- Comes in back and white

- SIM card provided Includes activation and \$35 credit applied
- Device programmed by us before despatch
- · Easily re-programmable
- Professional monitoring available
- Mobile phone live tracking/ monitoring APP available (Paid subscription)
- · Light weight
- Overall Battery Life: approx.
 5 years minimum
- Rechargeable battery & low battery reporting
- 2 3 days battery life when in Power Saving mode
- Waterproof to IP67

WANDERING ALERT



PROTRAC ID WANDERING ALERT SYSTEM

The Australian made system is complimentary to all nurse call systems in any size facility. Allowing up to 100 wrist tags to be worn at any one time the wandering alert system is discreet and versatile. The antenna which resembles a piece of Perspex can be mounted above the door and painted if desired otherwise placed inside the ceiling cavity. The Protrac wrist tags have been designed with longevity in mind. Using a low powered radio frequency. The wrist pendants will have a 3-5-year life span, unlike other wandering systems where wandering wrist tags are disposable after the first year and are more expensive and complex to program. The Protrac ID wandering system must be installed by a qualified installer – Contact us if you need assistance.



A-PROTRAC-WRIST TAG



PROTRAC ID WRIST TAGS/TRANSMITTERS

Protrac wrist tags

INCONTINENCE + SEIZURE DETECTION



CURA1 INCONTINENCE DETECTION BED SENSOR KIT

Used to detect and alert the carer when a loved one or a resident has been incontinent in bed..

KIT CONTAINS:

- 1 x incontinence bed pad 30 day use (including con-necting cord)
- · 1 x incontinence monitor
- · Replacement pads available



CURA1 INCONTINENCE DETECTION CHAIR SENSOR KIT

Used to detect and alert the carer when a loved one or a resident has been incontinent in a chair or wheelchair.

KIT CONTAINS:

- 1 x incontinence chair pad 30 day use (including con-necting cord)
- · 1 x incontinence monitor
- Replacement pads available



CURA1 SEIZUREALERT

Detect & raise an alarm as a seizure occurs. The SeizureAlert can detect seizure movement and also sound. The Cura1 SeizureAlert is easily installed and economical (no ongoing costs or consumables). With battery backup and minimal false alarms, it Connects easily with third party alarm devices and is able to be plugged into a nurse call system's bedside call point for nurse call connectivity.



With IP-Hardwired, ISM-Wireless, Bluetooth controlled RTLS, Presence, Duress, Dementia Care management, extensive reporting and accountability management.



WE ARE YOUR SOLUTION!

For a Nurse Call quotation or to discuss your requirements, please call the **Smart-Caller** sales team or visit the **Smart-Caller** website for more information.



www.smartcaller.com.au

(03) 9588 0833

sales@smartcaller.com.au